

charge domain global shutter

Pyramid Imaging

deep cooled down to -20° C

> compact design

appropriate for long exposure times

> high resolution 5120 x 5120 pixel





945 East 11th Avenue Tampa, FL 33605 Phone: (813) 984-0125 Contact: Sales@ pyramidimaging.com https://pyramidimaging.com



An Excelitas Technologies Brand


pco.edge 26

» sCMOS image sensor

type of sensor	global shutter scientific CMOS (sCMOS) monochrome
resolution (h x v)	5120 x 5120 active pixel
pixel size (h x v)	2.5 μm x 2.5 μm
sensor format / diagonal	12.8 mm x 12.8 mm / 18.1 mm
shutter mode	global/snapshot shutter (GS)
MTF	200 lp/mm (theoretical)
fullwell capacity	4500 e-
readout noise (typ.) ¹	2.3 med e ⁻ / 2.5 rms e ⁻
dynamic range (typ.)	66 dB
quantum efficiency	72 %
spectral range	320 nm 1000 nm
dark current (typ.)	0.09 e ⁻ /pixel/s @ -10 °C
DSNU	< 1 e ⁻ rms
PRNU	< 0.9 %
anti blooming factor ²	> 10 000
parasitic light sensitivity	1/10 000

» camera system

maximum frame rate @ full resolution	6 fps
exposure / shutter time	27 μs 60 s
dynamic range A/D	12 bit
A/D conversion factor	1.1 e ⁻ /DN
pixel data rate	187 Mpixel/s
binning horizontal	x1, x2, x4
binning vertical	x1, x2, x4
region of interest (ROI)	horizontal: steps of 8 pixel (min. 24) vertical: steps of 2 pixel (min. 8)
non linearity	< 0.6 %
cooling method	adjustable: from - 20 °C to + 20 °C peltier with forced air (fan) and water cooling calibration setpoint: - 10 °C
trigger input signals	frame trigger, sequence trigger, programmable input (SMA connectors)
trigger output signals	exposure, busy, programmable output (SMA connectors)
data interface	USB 3.1 Gen 1
time stamp	in image (1 µs resolution)

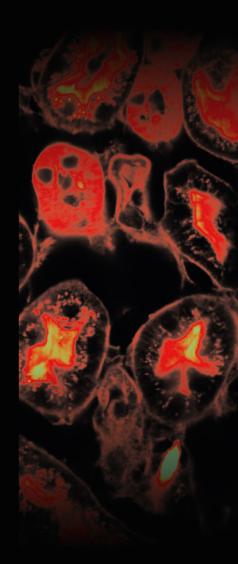

¹ The readout noise values are given as median (med) and root mean square (rms) values, due to the different noise models, which can be used for evaluation. All values are raw data without any filtering. ² Based on image sensor data sheet.



945 East 11th Avenue Tampa, FL 33605 Phone: (813) 984-0125 Contact: Sales@ pyramidimaging.com https://pyramidimaging.com

Pyramid Imaging


subject to changes without prior notice | lens is sold separately © Excelitas PCO GmbH, Kelheim | pco.edge 26 data sheet | v1.08




» general

power delivery	power over USB 3.1 Gen 1 and power connector (24 VDC +/- 10 %)
power consumption	22 W max. (typ. 10 W) over power connector and typ. 4.5 W over USB 3.1 Gen 1
weight	970 g
operating temperature	+ 10 °C + 40 °C
operating humidity range	10 % 80 % (non-condensing)
storage temperature range	- 10 °C + 60 °C
optical interface	F-mount, C-mount
maximum cable length	5 m
CE / FCC certified	yes

>> quantum efficiency

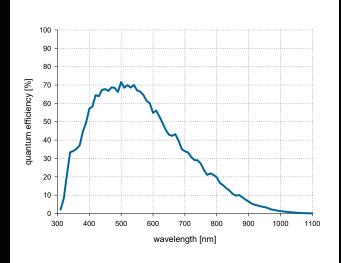

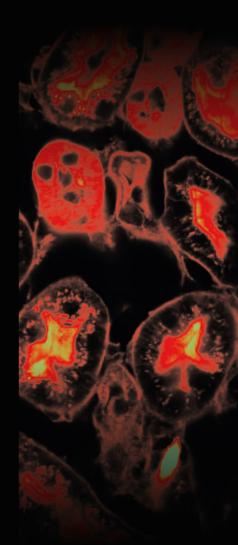

» frame rate table

5120 x 5120	6 fps
5120 x 1024	30 fps
5120 x 512	59 fps
5120 x 256	115 fps
5120 x 128	216 fps
1920 x 1080	29 fps
1600 x 1200	26 fps
1280 x 1024	30 fps
640 x 480	63 fps
320 x 240	122 fps




945 East 11th Avenue Tampa, FL 33605 Phone: (813) 984-0125 Contact: Sales@ pyramidimaging.com https://pyramidimaging.com



» benefits

anti-vibration technology with AlphaGEL® sealed electronics for dust & dirt protection optimized airflow with 3 seperate intakes spider-less mounting with only 6.18 mm to image plane

» dimensions

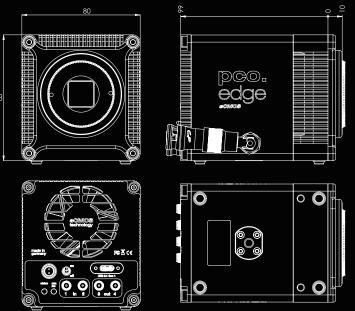

F-mount and C-mount lens adapter are changeable. All dimensions are given in millimeter.

>> camera view



Phone: (813) 984-0125

https://pyramidimaging.com







>> applications

brightfield microscopy | fluorescence microscopy | digital pathology | light sheet fluorescence microscopy (LSFM) | structured illumination microscopy (SIM) | high-speed bright field ratio imaging | high throughput screening | high content screening | biochip reading | spinning disk confocal microscopy | 3D metrology | industrial quality inspection

» software



With pco.camware you control all camera settings, the image acquisition and the storage of your image data. The pco.sdk is the complementary software development kit. It includes dynamic link libraries for user customization and integration on Windows-PC platforms. Drivers for popular third party software packages are also available for you.

All this items like pco.camware, pco.sdk and third party drivers, are free-to-download at <u>www.pco.de.</u>





contact

pco europe +49 9441 2005 50 info@pco.de pco.de

pco america +1 866 678 4566 info@pco-tech.com pco-tech.com



+65 6549 7054 info@pco-imaging.com pco-imaging.com

pco china

+86 512 67634643 info@pco.cn pco.cn



for application stories please visit our website

